

PRODUCT INFORMATION

Tag	C-Flag&Strep Tag
Target	HCAR2
Synonyms	GPR109A; HCA2; HM74a; HM74b; NIACR1; Puma-g; PUMAG
Description	Human HCAR2-Strep full length protein-synthetic nanodisc
Delivery	In Stock
Uniprot ID	Q8TDS4
Expression Host	HEK293
Protein Families	Druggable Genome, GPCR, Transmembrane
Protein Pathways	N/A
Molecular Weight	The human full length HCAR2-Strep protein has a MW of 41.9 kDa Lyophilized from nanodisc solubilization buffer (20 mM Tris-HCl, 150 mM NaCl, pH 8.0). Normally 5% - 8% trehalose is added as protectants before lyophilization. Please see Certificate of Analysis for specific instructions. Do not use solvents with a pH below 6.5 or those containing high concentrations of divalent metal ions (greater than 5 mM) in subsequent experiments.
Formulation & Reconstitution	Store at -20°C to -80°C for 12 months in lyophilized form. After reconstitution, if not intended for use within a month, aliquot and store at -80°C (Avoid repeated freezing and thawing). Lyophilized proteins are shipped at ambient temperature.
Storage & Shipping	Acts as a high affinity receptor for both nicotinic acid (also known as niacin) and (D)-beta-hydroxybutyrate and mediates increased adiponectin secretion and decreased lipolysis through G(i)-protein-mediated inhibition of adenylyl cyclase. This pharmacological effect requires nicotinic acid doses that are much higher than those provided by a normal diet. Mediates nicotinic acid-induced apoptosis in mature neutrophils. Receptor activation by nicotinic acid results in reduced cAMP levels which may affect activity of cAMP-dependent protein kinase A and phosphorylation of target proteins, leading to neutrophil apoptosis. The rank order of potency for the displacement of nicotinic acid binding is 5-methyl pyrazole-3-carboxylic acid = pyridine-3-acetic acid > acifran > 5-methyl nicotinic acid = acipimox >> nicotinuric acid = nicotinamide.
Background	
Usage	Research use only
Conjugate	Unconjugated



ELISA assay to evaluate HCAR2-Strep-Nanodisc 0.2 μ g Human HCAR2-Strep-Nanodisc per well

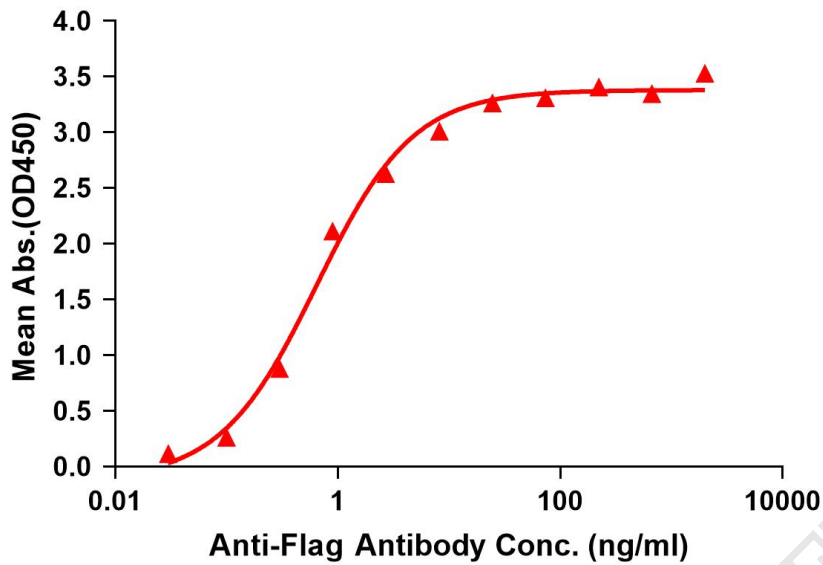


Figure 1. Elisa plates were pre-coated with C-Flag&Strep Tag HCAR2-Strep-Nanodisc (0.2 μ g/per well). Serial diluted anti-Flag monoclonal antibody solutions were added, washed, and incubated with secondary antibody before Elisa reading. From above data, the EC50 for anti-Flag monoclonal antibody binding with HCAR2-Strep-nanodisc is 0.649ng/ml.

kDa M R

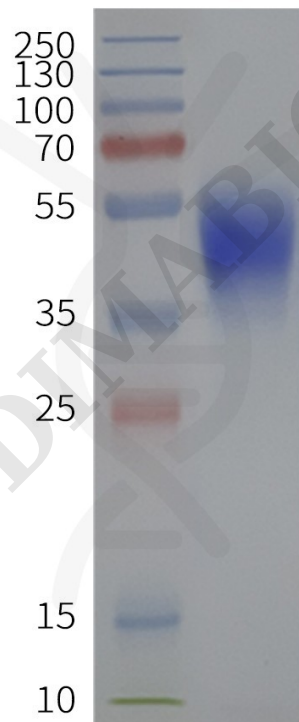


Figure 2. Human HCAR2-Strep-Nanodisc, C-Flag&Strep Tag on SDS-PAGE

